



Surgical Endoscopes


TRADITIONAL RIGID SCOPES

BR Surgical guarantees a German-engineered endoscope with premium image brightness and optical clarity that is designed and constructed for long lasting performance. Our service team is dedicated to keeping you in business with fast, fair, and cost-effective service solutions.

High quality construction and superior optics make this scope a long-lasting value. Excellent customer service and support make it a smart purchase.



Clearly a Better Scope

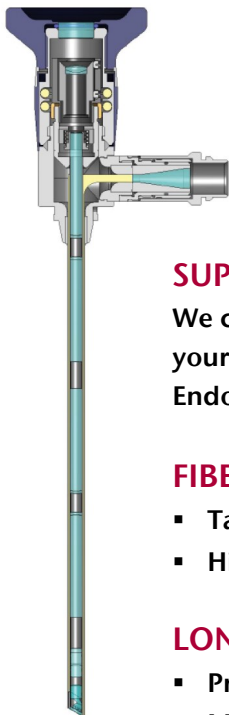

SUPERIOR VISUALIZATION & DURABILITY:

We can guarantee a brilliant high definition image because we control the construction of your endoscope from start to finish, from optical calculation to glass rod production. BR Endoscopes provide maximum protection against degrading influences of reprocessing.

FIBER OPTICS & ILLUMINATION:

- Tapered fiber cone = better light amplification
- High transmission wide angle illumination fibers = superior light transmission

LONG LASTING PERFORMANCE:

- Proximal & distal sapphire window = scratch resistant, no glass corrosion
- Medical grade stainless steel tubing = stronger, more durable material
- No coatings = eliminates chance of coatings coming off during reprocessing over time
- Laser welded joints = seals fiberoptic illumination bundles hermetically from the optical train, to avoid humidity and chemical penetration
- Laser welded joints = increase strength between scope body and sheath
- Hydrophilic gas-filled optical train = anti fogging device



All BR rigid endoscopes are autoclavable and come with a one year warranty against materials and workmanship

OVER 20,000 SURGICAL INSTRUMENTS & PROCEDURE SETS • ENDOSCOPIC VIDEO EQUIPMENT & RELATED MEDICAL DEVICES

BR SURGICAL • 3500 Beachwood Court • Suite 107 • Jacksonville, FL 32224 • customerservice@brsurgical.com
Office: 904.642.1366 • Toll Free: 888.642.1366 • Fax: 904.642.1368 • www.BRSurgical.com





Surgical Endoscopes


TRADITIONAL RIGID SCOPES

Rigid Endoscopes		Hysteroscopes (302mm length)	Cystoscopes (302mm length)	ENT/Arthroscopes (187/175/150mm length)
Ø 2.7mm	0°	-	-	BR946-2700-187
	30°	-	-	BR946-2730-187
	70°	-	-	BR946-2770-187
Ø 2.9mm	0°	BR970-2900-302	BR966-2900-302	-
	12°	BR970-2912-302	BR966-2912-302	-
	30°	BR970-2930-302	BR966-2930-302	-
Ø 3.0mm	0°	-	-	BR946-3000-150
	30°	-	-	BR946-3030-150
	45°	-	-	BR946-3045-150
	70°	-	-	BR946-3070-150
Ø 4.0mm	0°	BR970-4000-302	BR966-4000-302	BR946-4000-175
	12°	BR970-4012-302	BR966-4012-302	-
	30°	BR970-4030-302	BR966-4030-302	BR946-4030-175
	45°	-	-	BR946-4045-175
	70°	BR970-4070-302	BR966-4070-302	BR946-4070-175

Laparoscopes		0°	30°	45°
Ø 5.0mm	175mm length	BR953-5000-175	BR953-5030-175	-
	302mm length	BR953-5000-302	BR953-5030-302	-
	450mm length	BR953-5000-450	BR953-5030-450	BR953-5045-450
Ø 10.0mm	180mm length	BR953-1000-180	BR953-1030-180	-
	302mm length	BR953-1000-302	BR953-1030-302	BR953-1045-302
	450mm length	BR953-1000-450	BR953-1030-450	BR953-1045-450

Operating Laparoscopes

BR953-1055-277	Laparoscope, autoclavable, Ø 10mm, 0°, with 5.5mm working channel, 277mm
BR953-1060-270	Laparoscope, autoclavable, Ø 10mm, 0°, with 6mm working channel, 270mm

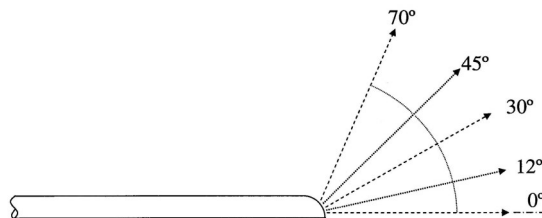

OVER 20,000 SURGICAL INSTRUMENTS & PROCEDURE SETS • ENDOSCOPIC VIDEO EQUIPMENT & RELATED MEDICAL DEVICES

BR SURGICAL • 3500 Beachwood Court • Suite 107 • Jacksonville, FL 32224 • customerservice@brsurgical.com
 Office: 904.642.1366 • Toll Free: 888.642.1366 • Fax: 904.642.1368 • www.BRSurgical.com

